



参数规格SPECIFICATION

DISPOSABLE BIOPSY FORCEPS 一次性活体取样钳

规格型号 Item No.	钳子杯直径 O.D.(mm)	有效长度 Working Length (mm)	钳口形式 Jaw Shape	是否带针 Spike(N)	是否包塑 Coating(B)	最小可通过钳道 Endoscope Channel (mm)
PTY1-23(16/18/20/23) T	2.3	1600/1800/2000/2300	Oval平口(T)	No	No	2.8
PTY1-23(16/18/20/23) TN	2.3	1600/1800/2000/2300	Oval平口(T)	Yes	No	2.8
PTY1-23(16/18/20/23) TB	2.3	1600/1800/2000/2300	Oval平口(T)	No	Yes	2.8
PTY1-23(16/18/20/23) TNB	2.3	1600/1800/2000/2300	Oval平口(T)	Yes	Yes	2.8
PTY1-23(16/18/20/23)U	2.3	1600/1800/2000/2300	Alligator鳄鱼口(U)	No	No	2.8
PTY1-23(16/18/20/23)UN	2.3	1600/1800/2000/2300	Alligator鳄鱼口(U)	Yes	No	2.8
PTY1-23(16/18/20/23)UB	2.3	1600/1800/2000/2300	Alligator鳄鱼口(U)	No	Yes	2.8
PTY1-23(16/18/20/23)UNB	2.3	1600/1800/2000/2300	Alligator鳄鱼口(U)	Yes	Yes	2.8
PTY1-23(16/18/20/23)T-O	2.3	1600/1800/2000/2300	Oval平口(T), Distal fine对接 (O)	No	No	2.8
PTY1-23(16/18/20/23)TN-O	2.3	1600/1800/2000/2300	Oval平口(T), Distal fine对接 (O)	Yes	No	2.8
PTY1-23(16/18/20/23)TB-M	2.3	1600/1800/2000/2300	Oval平口(T), Distal fine对接 (O)	No	Yes	2.8
PTY1-23(16/18/20/23)TNB-M	2.3	1600/1800/2000/2300	Oval平口(T), Distal fine对接 (O)	Yes	Yes	2.8
PTY1-23(16/18/20/23)TB-D	2.3	1600/1800/2000/2300	Oval平口(T), Unilateral turn 单侧翻 (D)	No	Yes	2.8
PTY1-23(16/18/20/23)TB-S	2.3	1600/1800/2000/2300	Oval平口(T), Bilateral turn 双侧翻 (S)	No	Yes	2.8
PTY2-18(07/10/12/16/18/23)T	1.8	700/1000/1200/1600/1800/2300	Oval平口(T)	No	No	2.0
PTY2-18(07/10/12/16/18/23)TB	1.8	700/1000/1200/1600/1800/2300	Oval平口(T)	No	Yes	2.0
PTY2-18(07/10/12/16/18/23)U	1.8	700/1000/1200/1600/1800/2300	Alligator鳄鱼口(U)	No	No	2.0
PTY2-18(07/10/12/16/18/23)UB	1.8	700/1000/1200/1600/1800/2300	Alligator鳄鱼口(U)	No	No	2.0

Reusable Biopsy Forceps 重复性活检钳

Item No.	O.D.(mm)	Working Length(mm)	Jaw Shape	Spike	Coating	Endoscope Channel(mm)
PTRF-18(16/18/23)T	1.8	1600/1800/2300	Oval平口(T)	No	No	2.0
PTRF-18(16/18/23)U	1.8	1600/1800/2300	Alligator鳄鱼口(U)	No	No	2.0
PTRF-23(16/18/23)T	2.3	1600/1800/2300	Oval平口(T)	No	No	2.8
PTRF-23(16/18/23)U	2.3	1600/1800/2300	Alligator鳄鱼口(U)	No	No	2.8

产品特点Product Features

用于内窥镜下消化道、呼吸道活体组织取样

This device is used for tissue sampling in the digestive tract or respiratory tract through endoscope.

The alligator teeth applies for hard tissue sampling and takes an easily tumor tissue.

With needle it can help accurate sampling.

Coating outside effectively decreases the friction to the endoscope channel.

Moving cups are easy for accurate targeting, no sliding, depressing lesions.



诸暨市鹏天医疗器械有限公司
ZHUJI PENGtian MEDICAL INSTRUMENT CO.,LTD.

电话Tel:86-575-87128772 手机Mob:86-18657572931 (INT Market国外销售) 86-13777302827(China 国内销售) 传真Fax:86-87128771 网址:www.pengtianmed.com

邮箱Email:allen@pengtianmed.com 地址:中国浙江诸暨暨阳经济开发区津津路8号 Add:No.8 Jinjin Road,Jiyang Economic Development District311800,Zhuji,Zhejiang,China



参数规格SPECIFICATION

DISPOSABLE POLYPECTOMY SNARES一次性使用电圈套器

Item No. 规格型号	O.D. 外管直径 (mm)	Working Length 有效长度(mm)	Loop Diameter 圈径(mm)	Loop Type 圈套器形状	Endoscope Channel 最小可通过钳道(mm)
PTS-18(12/18)-(15/20/25)-(A/B)	1.8	1200/1800	15/20/25	Oval椭圆形(A) Hexagonal六角形(B)	2.0
PTS-23(16/18/23)-(10/15/20/25/30/35/40)-(A/B)	2.3	1600/1800/2300	10/15/20/25/30/35/40	Oval椭圆形(A) Hexagonal六角形(B)	2.8
PTS-23(16/18/23)-(10/15/20/25)-C	2.3	1600/1800/2300	10/15/20/25	Crescent半月形(C)	2.8

DISPOSABLE MONOFILAMENT SNARES一次性使用单股丝圈套器

Item No. 规格型号	O.D. 外管直径(mm)	Working Length 有效长度(mm)	Loop Diameter 圈径(mm)	Loop Type 圈套器形状	Endoscope Channel 最小可通过钳道(mm)
PTS-23(16/18/23)-(10/15/20/25/30)-(A/B)-D	2.3	1600/1800/2300	10/15/20/25/30	Oval椭圆形(A) Hexagonal六角形(B)	2.8
PTS-23(16/18/23)-(10/15/20/25)-C-D	2.3	1600/1800/2300	10/15/20/25	Crescent半月形(C)	2.8

DISPOSABLE COLD SNARES一次性使用冷圈套器

Item No. 规格型号	O.D. 外管直径(mm)	Working Length 有效长度(mm)	Loop Diameter 圈径(mm)	Loop Type 圈套器形状	Endoscope Channel 最小可通过钳道(mm)
PTS-23(16/18/23)-(6/10/15)-(A/B/E)-L	2.3	1600/1800/2300	6/10/15	Oval椭圆形(A) Hexagonal六角形(B) Half Hexagonal(E)	2.8

DISPOSABLE ROTATED SYNCHRONOUSLY SNARES同步旋转圈套器

Item No. 规格型号	O.D. 外管直径(mm)	Working Length 有效长度(mm)	Loop Diameter 圈径(mm)	Loop Type 圈套器形状	Endoscope Channel 最小可通过钳道(mm)
PTS-23(16/18/23)-(10/15/20/25/30/35/40)-A-T	2.3	1600/1800/2300	10/15/20/25/30/35/40	Oval椭圆形(A)	2.8
PTS-23(16/18/23)-(10/15/20/25/30/35)-B-T	2.3	1600/1800/2300	10/15/20/25/30/35	Hexagonal六角形(B)	2.8

适用于直产品特点Product features

电圈套器与高频电流配合使用，用于经内窥镜对消化道、呼吸道的息肉的切除

This polypectomy snare is used with high frequency to resect small polyp under endoscope.

为适应不同息肉形状的需要，有椭圆形、六角形、半六角形和半月形不同套圈形状

Oval, hexagonal, half-hexagonal and crescent type of snare are supplied to meet the needs of different shapes of polyp.

连接牢固、可旋转操作，操控性好，可实现单手操作

Connect firmly, available for rotation operation, and well controllability, can be operated by single hand.

单股丝圈套器由镍钛合金丝组成，能够保持很好的套圈形状。锋利的套圈，可以完整地切除息肉，留下很小的创伤面积

It's made of nickel-titanium alloy wire. The shape of loop can be maintained in a good status. With sharp loop. The polyps can cut completely while with rare wound area

特别设计的薄套圈适用于切割直径≤10mm的息肉，并切能够完整的切除

Specially-designed thin loop wire cold snare allows for a clean, guillotine type cut and it is suitable for polyps ≤10mm in size

同步旋转圈套器的套圈是由高强度医用编织钢丝制成，具有精确快速的切割性能，同时也具有优异的电凝性能。可同步旋转特性非常便于操作



The rotatable synchronously snare is made of high-strength medical braided wire, offering a precise & quick cutting properties, as well as an outstanding electric coagulation. The rotatable synchronously characteristic is very convenient for operation



诸暨市鹏天医疗器械有限公司
ZHUJI PENG TIAN MEDICAL INSTRUMENT CO.,LTD.

电话Tel:86-575-87128772 手机Mob:86-18657572931 (INT Market国外销售) 86-13777302827(China 国内销售) 传真Fax:86-87128771 网址:www.pengtianmed.com

邮箱Email:allen@pengtianmed.com 地址:中国浙江诸暨暨阳经济开发区津津路8号 Add:No.8 Jinjin Road,Jiyang Economic Development District311800,Zhuji,Zhejiang,China

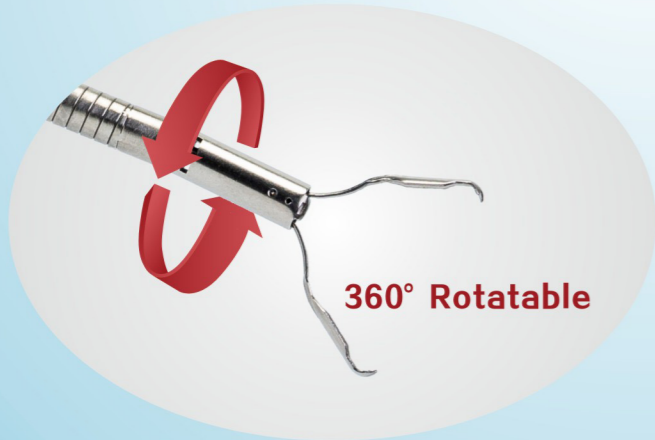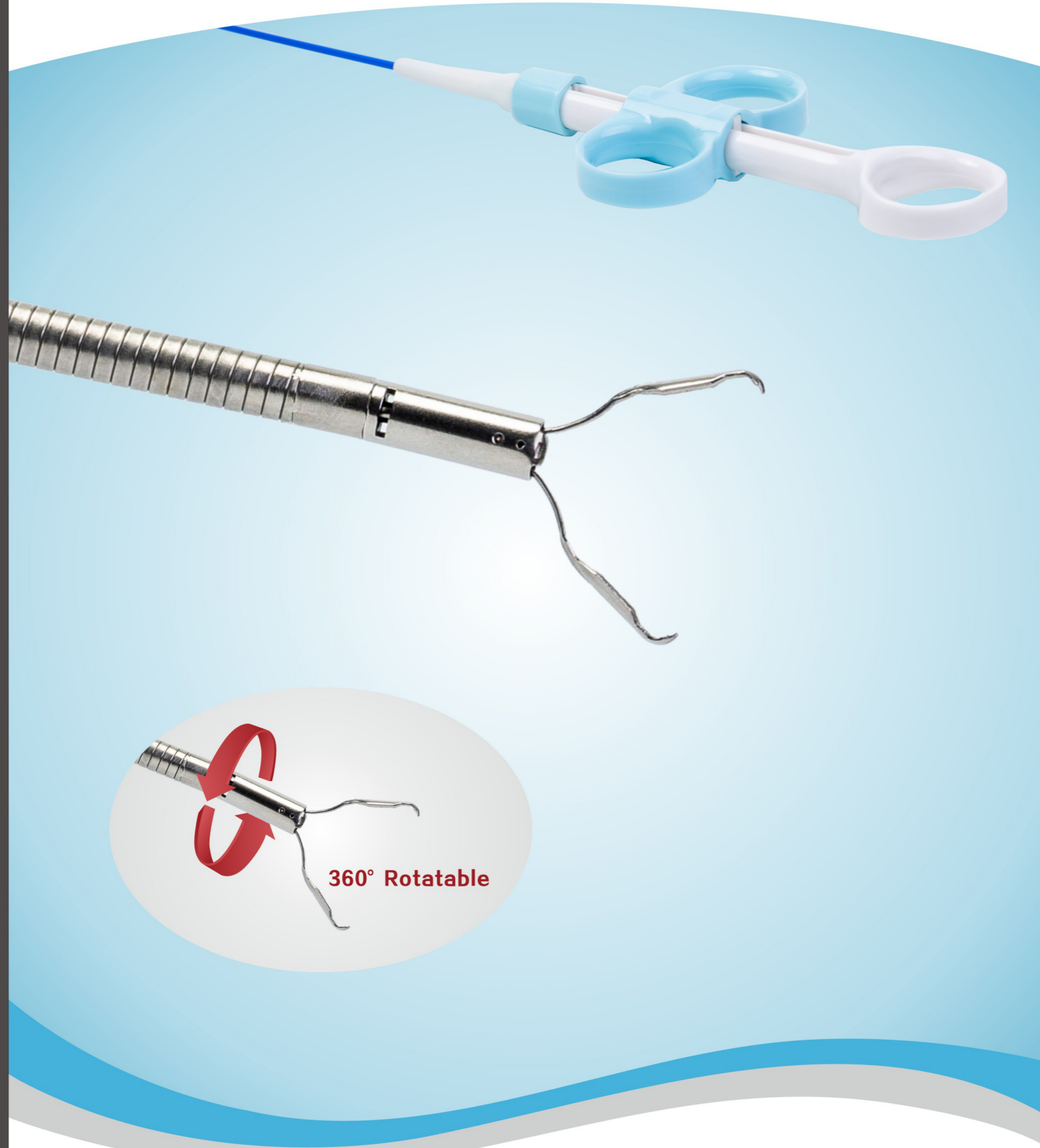

参数规格SPECIFICATION



DISPOSABLE HEMOCLIPS 夹子装置

Item No. 规格型号	Opening Size 开口尺寸 (mm)	O.D. 外管直径 (mm)	Working Length 工作长度(mm)	Coating 是否包塑
PTZJ-(C/D/E/F)-26-(160/170/180/190/195/200/210/ 220/230/270)	C≥8 D≥10 E≥12 F≥14	2.6	1600/1700/1800/1900/1950/ 2000/2100/2200/2300/2700	No
PTZJ-(C/D/E/F)-26-(160/170/180/190/195/200/210/ 220/230/270)-B	C≥8 D≥10 E≥12 F≥14	2.6	1600/1700/1800/1900/1950/ 2000/2100/2200/2300/2700	Yes

产品特点Product features

用于在内窥镜下夹合消化道内软组织，消化道止血、标记

This device is used for clamping soft tissue, gastrointestinal bleeding, marking under endoscope.

可重复开闭，真正实现重复预夹功能；

The clips can be opened and closed repeatedly, which truly realize the function of pre-clamping.

可旋转功能，有助于医生以最合适的角度夹持组织。

This clips can be rotated to help the doctor hold the tissue at the most appropriate angle.

Different sizes suit different needs of adults and children.



诸暨市鹏天医疗器械有限公司
ZHUJI PENG TIAN MEDICAL INSTRUMENT CO.,LTD.

电话Tel:86-575-87128772 手机Mob:86-18657572931 (INT Market国外销售) 86-13777302827 (China 国内销售) 传真Fax:86-87128771 网址:www.pengtianmed.com
邮箱Email:allen@pengtianmed.com 地址:中国浙江诸暨暨阳经济开发区津津路8号 Add:No.8 Jinjin Road,Jiyang Economic Development District311800,Zhuji,Zhejiang,China



DISPOSABLE ENDOSCOPE INJECTION NEEDLES 一次性内镜用注射针

Item No. 规格型号	O.D. 外管直径(mm)	Working Length 工作长度(mm)	Needle Diameter 针径 (mm)	Needle Length 针长 (mm)
PTI-18(16/18/23)-(22/23/25) (04/05/06)-C	Φ1.8	1600/1800/2300	22G(0.7) 23G(0.6) 25G(0.5)	4/5/6
PTI-23(16/18/23)-(22/23/25) (04/05/06)-(A/B/C)	Φ2.3	1600/1800/2300	22G(0.7) 23G(0.6) 25G(0.5)	4/5/6


产品特点 Product features

- 手柄单手操作，方便医生临床操作
One-handed performance, which is convenient for doctors to operate clinically.
- 管身透明，回血清晰，安全保证
Transparent tube body, clear blood return, safety guarantee.
- 针尖锐利，进针顺滑，易于穿刺
Sharp needle tip, smooth needle insertion, easy to puncture
- 管身传导性好，弯曲盘旋后仍可顺利出针
The tube body has good conductivity, and the needle can still be smoothly withdrawn after bending and circling



诸暨市鹏天医疗器械有限公司
ZHUJI PENG TIAN MEDICAL INSTRUMENT CO.,LTD.

电话Tel:86-575-87128772 手机Mob:86-18657572931 (INT Market国外销售) 86-13777302827(China 国内销售) 传真Fax:86-87128771 网址:www.pengtianmed.com
邮箱Email:allen@pengtianmed.com 地址:中国浙江诸暨暨阳经济开发区津津路8号 Add:No.8 Jinjin Road,Jiyang Economic Development District311800,Zhuji,Zhejiang,China

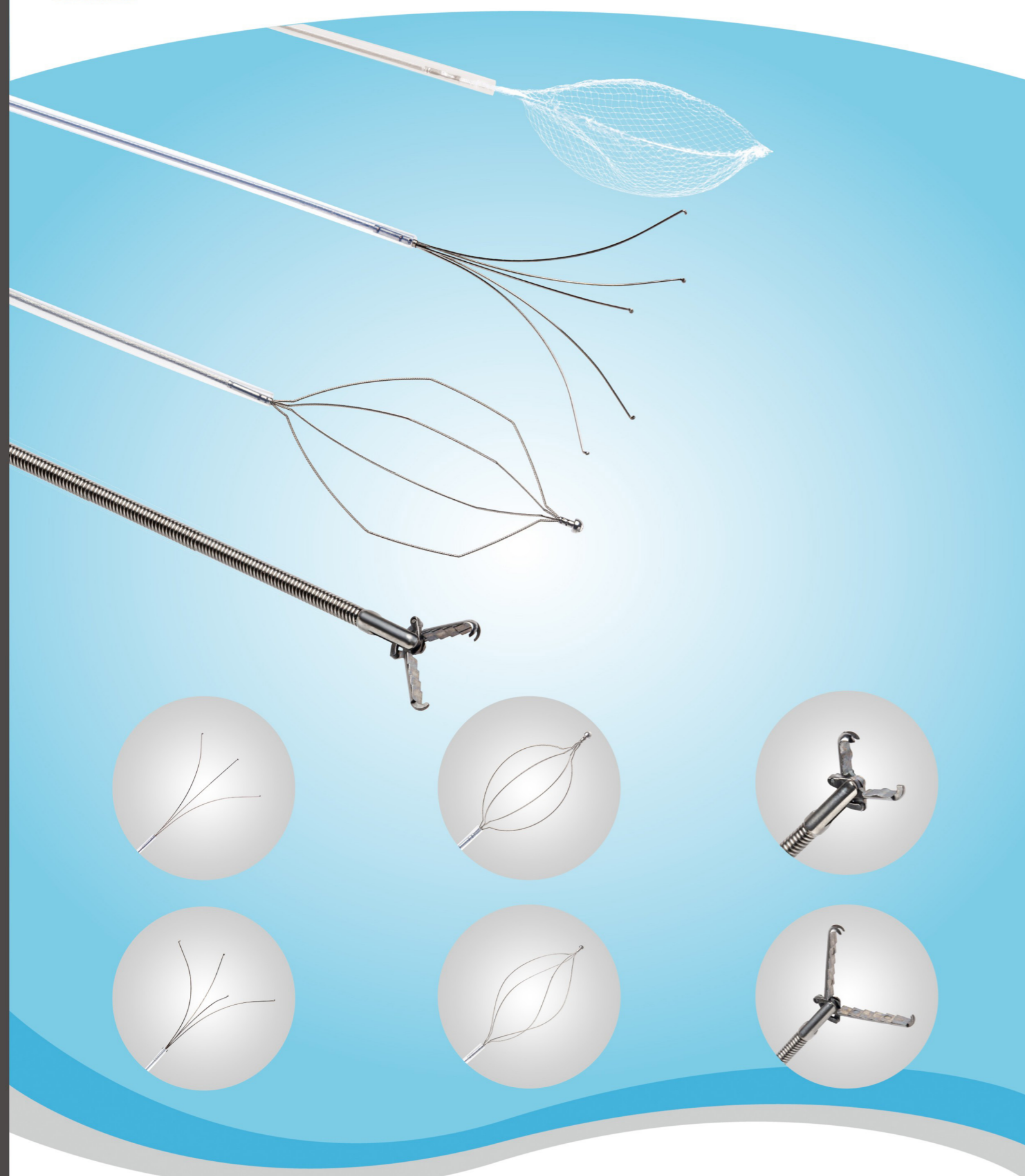

参数规格SPECIFICATION

DISPOSABLE GRASPING FORCEPS异物钳

Item No. 规格型号	O.D. 外管直径(mm)	Working Length 有效长度(mm)	Coating 包塑	Prong/Jaw Shape 爪脚数量/钳口形式
PTG-2316- (Z3/Z4/Z5)	2.3	1600/1800/2300	/	3/4/5 prongs
PTG-18(16/18/23)C	1.8	1600/1800/2300	No	Rat Tooth With Alligator Jaws鼠齿鳄鱼齿型
PTG-23(16/18/23)S	2.3	1600/1800/2300	No	Rat Tooth Type鼠齿型
PTG-23(16/18/23)SB	2.3	1600/1800/2300	Yes	Rat Tooth Type鼠齿型
PTG-23(16/18/23)E	2.3	1600/1800/2300	No	Alligator Type鳄鱼齿型
PTG-23(16/18/23)EB	2.3	1600/1800/2300	Yes	Alligator Type鳄鱼齿型
PTG-23(16/18/23)C	2.3	1600/1800/2300	No	Rat Tooth With Alligator Jaws鼠齿鳄鱼齿型
PTG-23(16/18/23)CB	2.3	1600/1800/2300	Yes	Rat Tooth With Alligator Jaws鼠齿鳄鱼齿型
PTG-23(16/18/23)T	2.3	1600/1800/2300	No	Pelican Type 塘鹅嘴型
PTG-23(16/18/23)TB	2.3	1600/1800/2300	Yes	Pelican Type 塘鹅嘴型

REUSABLE GRASPING FORCEPS重复性异物钳

Item No. 规格型号	O.D. 外管直径(mm)	Working Length 有效长度(mm)	Prong/Jaw Shape 爪脚数量/钳口形式
PTRG-18(16/18/23)C	1.8	1600/1800/2300	Rat Tooth With Alligator Jaws鼠齿鳄鱼齿型
PTRG-23(16/18/23)S	2.3	1600/1800/2300	Rat Tooth Type鼠齿型
PTRG-23(16/18/23)E	2.3	1600/1800/2300	Alligator Type鳄鱼齿型
PTRG-23(16/18/23)C	2.3	1600/1800/2300	Rat Tooth With Alligator Jaws鼠齿鳄鱼齿型
PTRG-23(16/18/23)T	2.3	1600/1800/2300	Pelican Type 塘鹅嘴型

用于经内窥镜抓取切下的息肉、病变组织、团状异物
This device is used for grasping the cut polyps, lesion tissue or mass foreign body.

RETRIEVAL NET网兜

Item No. 规格型号	O.D. 外管直径(mm)	Working Length 有效长度(mm)	Opening Diameter 开幅直径(mm)
PTG-C-(20/25/30/35/40)-26(160/180/230)	2.6	1600/1800/2300	20/25/30/35/40

带网设计, 可准确牢固地抓取异物
This device can accurately grasp the foreign body firmly with a net design.

DISPOSABLE STONE EXTRACTION BASKETS一次性使用取石网篮

Item No. 规格型号	O.D. 外管直径(mm)	Working Length 工作长度(mm)	Basket Width 网篮开幅(mm)	Basket Shape 网篮类型	Endoscope Channel 最小可通过钳道(mm)
PTQ-18(07/10/18)-(20/25/30)-(A/B/C)	1.8	700/1000/1800	20/25/30	Diamond 钻石形(A) Oval 椭圆形(B) Spiral 螺旋形(C)	2.0
PTQ-23(16/18/23)-(20/25/30)-(A/B/C)	2.3	1600/1800/2300	20/25/30	Diamond 钻石形(A) Oval 椭圆形(B) Spiral 螺旋形(C)	2.8

用于抓取内窥镜下结石、息肉切除后的息肉组织或异物
This device is used to grasp the stones, polyps after removal of polyps or foreign bodies under endoscope.
旋转手柄可以实现网篮的360°旋转, 有效抓取小结石和粉碎性结石
Rotatable handle can make the basket a 360° rotation, which effectively captures small stones and crushed stones.
不锈钢网篮头拉紧定型能力好, 套取后不易脱落
Medical-use stainless steel tightens the basket well, which makes it uneasy to fall off after taking.





参数规格 SPECIFICATION

DISPOSABLE HOT BIOPSY 一次性使用热活检钳

Item No. 规格型号	O.D. 外管直径(mm)	Working Length 有效长度(mm)	Endoscope Channel 最小可通过钳道(mm)
PTD-J1/J2/J3/J4-23(07/10/12/16/18/23)	2.3	700/1000/1200/1600/1800/2300	2.8

与高频电流配合使用，用于内窥镜下消化道、呼吸道的切割和凝血
This device is used for polyp resection and blood coagulation in the digestive tract or respiratory tract through endoscope with high frequency.
快速电凝出血部位
Coagulate bleeding site rapidly
旋转功能，快速止血
Rotating function and Stop bleeding fast
安全抓住出血点
To securely grasp bleeding points.
外表涂层起绝缘作用及保护内窥镜不受损伤
Insulated coating outside effectively protects and decreases the friction to the endoscope channel.

BITE BLOCKS 口垫

Item No. 规格型号	Characteristics 特征	Item No. 规格型号	Characteristics 特征
PTB-YO-01	For adult, without belt	PTB-YO-05	For adult,with oxygen pipe, without belt
PTB-YO-02	For adult, with belt	PTB-YO-06	For pediatric,without belt
PTB-YO-03	For pediatric,with belt	PTB-YO-07	For adult,with belt and tongue
PTB-YO-04	For adult,with oxygen pipe and belt		

表面光滑，有效保护患者口腔和内窥镜钳道不受损伤
Smooth surface can effectively protect oral cavity of patients and avoid damage to endoscope channel.
不同规格适合大人小孩的不同需求
Different sizes suit different needs of adults and children.

BIOPSY VALVE 活检阀

Item No. 规格型号	Characteristics 特征
PT-HJF-01	For Olympus
PT-HJF-02	For Pentax

使内窥镜附件顺畅的进入钳道，并保护钳道，减少内窥镜损伤
This device is used to make endoscopic accessory smoothly enter the clamp and avoid injury of endoscope channel.



诸暨市鹏天医疗器械有限公司
ZHUJI PENG TIAN MEDICAL INSTRUMENT CO.,LTD.

电话Tel:86-575-87128772 手机Mob:86-18657572931 (INT Market国外销售) 86-13777302827(China 国内销售) 传真Fax:86-87128771 网址:www.pengtianmed.com
邮箱Email:allen@pengtianmed.com 地址:中国浙江诸暨暨阳经济开发区津津路8号 Add:No.8 Jinjin Road,Jiyang Economic Development District311800,Zhuji,Zhejiang,China



参数规格SPECIFICATION



DISPOSABLE ENDOSCOPE ULTRASOUND SUCTION BIOPSY NEEDLES
一次性内窥镜超声吸引活检针

Item No. 规格型号	Needle Diameter		O.D. 外管直径(mm)	Effective Length 有效长度 (mm)	Adjustable Length of Outer Tube 外管调节长度(mm)	Needle Adjustment Length 针管调节长度(mm)
	Gauge	针管直径(mm)				
PTN-A-S-10/16/1400/30/80	19G	1	Φ1.6	1400	30	80
PTN-A-S-08/16/1400/30/80	21G	0.8	Φ1.6	1400	30	80
PTN-A-S-07/16/1400/30/80	22G	0.7	Φ1.6	1400	30	80

产品特点Product features

与超声内窥镜配套使用，用于消化道进行超声引导下细针穿刺活检和组织取样。
Used in conjunction with an ultrasound endoscope for ultrasound-guided fine-needle aspiration biopsy and tissue sampling of the digestive tract.

独特的针管设计，穿刺性能好，可高效获取优质的细胞学/组织学样本。
Unique needle design, good puncture performance, can efficiently obtain high-quality cytology/histology samples.

外管可进行一定长度调节以匹配不同超声内窥镜。
The outer tube can be adjusted to a certain length to match different ultrasound endoscopes.

手柄防滑设计，管身顺滑性好，过弯性顺畅，灵活插入。
Anti-slip handle design, good smoothness, smooth bending, flexible insertion.



诸暨市鹏天医疗器械有限公司
ZHUJI PENG TIAN MEDICAL INSTRUMENT CO.,LTD.

电话Tel:86-575-87128772 手机Mob:86-18657572931 (INT Market国外销售) 86-13777302827(China 国内销售) 传真Fax:86-87128771 网址www.pengtianmed.com
邮箱Email:allen@pengtianmed.com 地址:中国浙江诸暨暨阳经济开发区津津路8号 Add:No.8 Jinjin Road,Jiyang Economic Development District311800,Zhuji,Zhejiang,China



参数规格(国内型号) SPECIFICATION (DOMESTIC MARKET)

DISPOSABLE DOUBLE ENDED CLEANING BRUSHES 一次性清洗刷

Item No. 规格型号	O.D. 外管直径(mm)	Working Length 有效长度(mm)	Brush Diameter 刷头直径(mm)
PTCB-A-03-16/18/23	1.8	1600/1800/2300	3
PTCB-A-06-16/18/23		1600/1800/2300	6
PTCB-A-06/11-16/18/23		1600/1800/2300	6/11

REUSABLE CLEANING BRUSHES 重复性清洗刷

Item No. 规格型号	O.D. 外管直径(mm)	Working Length 有效长度(mm)	Brush Diameter 刷头直径(mm)
PTCB-C-06-12/16/18/23-D	1.8	1200/1600/1800/2300	6
PTCB-C-06-12/16/18/23-S	1.8	1200/1600/1800/2300	6

用于清洗内窥镜钳道

Cleaning brushes are used to clean the endoscope channel.

一次性双头清洗刷使用方便, 高效清洗

Disposable double ended cleaning brushes are used conveniently and effectively.

一次性、重复性可供不同选择

Disposable type and reusable type are available for different choices.

DISPOSABLE CYTOLOGY BRUSHES 一次性使用细胞刷

Item No. 规格	O.D. 外管直径 (mm)	Working Length 工作长度(mm)	Brush Head Model 刷头型号	Brush Diameter 刷头直径(mm)	Bristle Length 刷长(mm)	Endoscope Channel 最小可通过钳道(mm)
PTX-18(10/12/16/18)-I-(20/30/40/50)10	1.8	1000/1200/1600/1800	Straight 直型刷头(I)	2/3/4/5	10	2.0
PTX-23(16/18/23)-I-(20/30/40/50)10	2.3	1600/1800/2300	Straight 直型刷头(I)	2/3/4/5	10	2.8

用于内窥镜下消化道呼吸道组织取样

The device is used for endoscopic tissue sampling of the digestive tract or respiratory tract.

球形光滑头端, 有效降低对腔壁的刺激和损伤

Spherical smooth head effectively reduces the stimulation and damage of the cavity wall.

直形细胞刷使得容易进入消化道呼吸道深处

Straight type makes it easy to enter the depth of the digestive tract or respiratory tract.

PORT HOLE CLEANING BRUSHES 水汽道刷

Item No. 规格型号	Brush Diameter(mm) 刷头直径
PTCB-B	5/11

用于内镜水汽按钮管道清洗。

The device is used for cleaning endoscopic water gas button pipe.



诸暨市鹏天医疗器械有限公司
ZHUJI PENG TIAN MEDICAL INSTRUMENT CO.,LTD.

电话Tel:86-575-87128772 手机Mob:86-18657572931 (INT Market国外销售) 86-13777302827(China 国内销售) 传真Fax:86-87128771 网址:www.pengtianmed.com

邮箱Email:allen@pengtianmed.com 地址:中国浙江诸暨暨阳经济开发区津津路8号 Add:No.8 Jinjin Road,Jiyang Economic Development District311800,Zhuji,Zhejiang,China